


위험한 EtBr은 이제 그만!!

Non-Mutagenic, Non-cytotoxic 형광다이


TopGreen & TopRed

Features


- Safety: Nonmutagenic and noncytotoxic
- Easy disposal: Safe to dispose in the drain
- Compatibility: Spectrally compatible with existing instruments
- Sensitivity: Higher signal but lower background
- Stability: can be stored at RT and microwavable

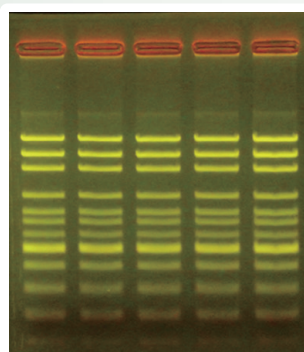


(1) TGN-001P

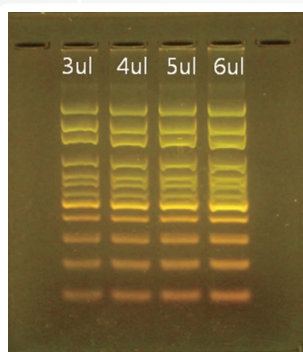

(2) TGN-002P

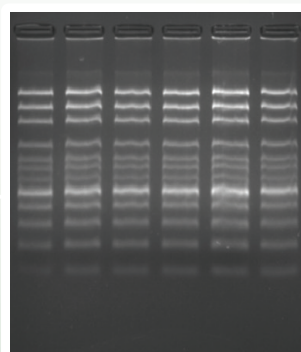

(3) TRN-003

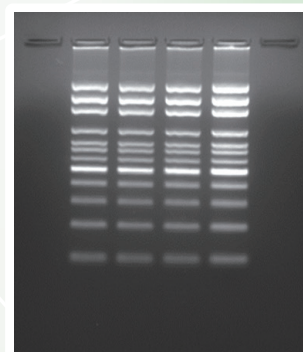

(4) TRN-004

- (1) : 2% agarose에 DNA sample 4ul에 TGN-001P 1ul를 sample staining하여 **LED** gel doc.으로 촬영
- (2) : 2% agarose 50ml에 TGN-002P 5ul, DNA loading 양 3,4,5,6ul를 전기영동 한 후 **LED** gel doc.으로 촬영
- (3) : 2% agarose에 DNA sample 4ul에 TRN-003 1ul를 sample staining하여 **UV** gel doc.으로 촬영
- (4) : 2% agarose 50ml에 TRN-004 5ul, DNA loading 양 3ul,5ul를 전기영동 한 후 **UV** gel doc.으로 촬영

Cat. No.	Name	Size	Price (별도가)
TGN-001P	TopGreen Nucleic acid 6X loading dye(LED)	1ml	50,000원
TGN-002P	TopGreen Nucleic acid Gel Stain(LED)	500ul	150,000원
TRN-003	TopRed Nucleic acid 6X loading dye (UV)	1ml	50,000원
TRN-004	TopRed Nucleic acid Gel Stain(UV)	500ul	150,000원

GENOMIC BASE
FOR YOUR COMFORT

(주)제노믹베이스

제품 문의 및 판매처

경기도 남양주시 다산동 6144-1
다산 진건 블루웨이브 지식 산업센터 2차 F315호
Tel: (02) 2215-4925 Fax: (02) 2215-4929
E-mail: info@genomicbase.com

지역 대리점

서울: 경기특약점: 와이젠바이오
강원: 세진사이언스
충북: (주)유티젠
충남대전: 하이바이오
전북: 제노셀

02) 858-7905
033) 764-4502
043) 266-5842
042) 825-7209
070) 4759-4197

경북대구: 지오리서치
전남: (주)동아하이텍
경남: 바이오드림(주)
부산: 대경바이오텍
제주: 비제이사이언스

053) 710-4790
062) 531-2701
055) 762-8059
010-5425-1780
064) 722-8711